


CE

# Aegistab OP850

## Electronic Operating Table

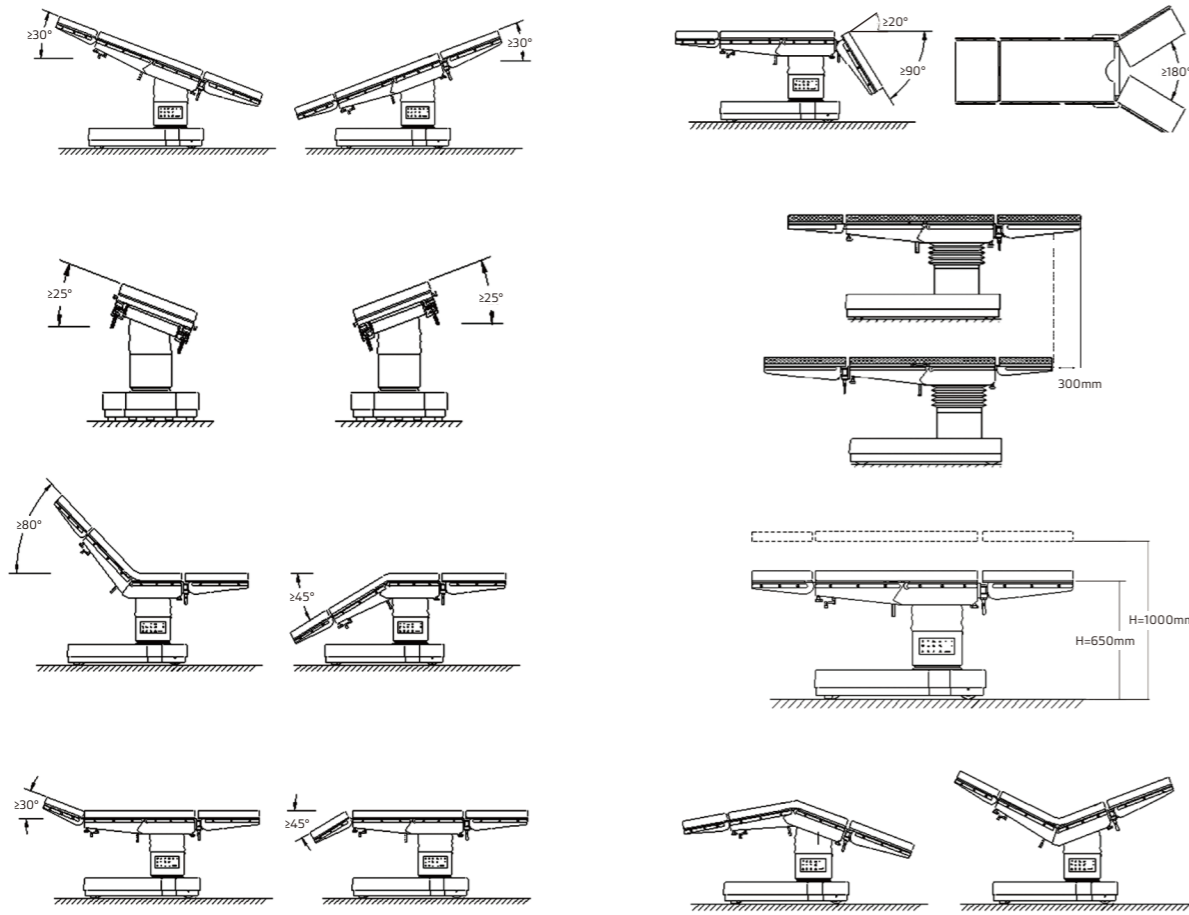

### Technical Specifications

Model	OP850 Standard	OP850 One Key Reset	OP850 One Key Reset, Shifting
Table Length (mm)	2080	2080	2080
Table Width (without side rails) (mm)	520	520	520
Table Height (without pads) (mm)	650 to 1000	650 to 1000	650 to 1000
Trendelenburg/Reverse	30°/30°	30°/30°	30°/30°
Lateral Tilt (left/right)	25°/25°	25°/25°	25°/25°
Back Plate Up/Down	80°/45°	85°/45°	85°/45°
Head Plate Up/Down (detachable)	30°/45°	30°/45°	30°/45°
Leg Plate Up/Down (detachable)	20°/90°	20°/90°	20°/90°
Leg Plate Split (detachable)	180°	180°	180°
Flex/Reflex	225° / 110°	225° / 110°	225° / 110°
Kidney-bridge Elevator Up (mm)	128	128	128
Longitudinal Shift (mm)	N/A	N/A	300
Weight Capacity (kg)	350	350	350
Exchangeable Head and Leg Plates	Yes	Yes	Yes
One Key Reset	N/A	Yes	Yes




Corporate Headquarter  
**HEYER Medical AG**  
 Carl-Heyer-Str. 1/3, 56130 Bad Embs, Germany  
 Tel: +49 2603 791 0  
 Fax: +49 2603 791 200  
 info@heyerm.com  
 www.heyerm.com

Manufacturer  
**HEYER Medical Co., Ltd.**  
 No.10, Chaobai Street, Yanjiao Development Zone  
 065201, Sanhe City, Hebei, China  
 Tel: +86 10 5841 1198  
 Fax: +86 10 6371 8989  
 marketing@heyerm.com



A high quality electronic-hydraulic operating table with a maximum load of 350kg, with versatile surgical configurations meeting the requirements of different surgeries.

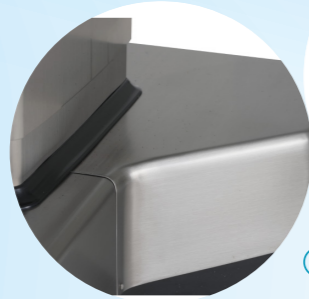

① Premium Quality Material

- All exposed parts are made of HI-Q 304 stainless steel
- High corrosion resistance, easy to clean and disinfection




② Built-in Elevator Kidney-bridge

- Built-in X-ray radiolucent kidney-bridge supports kidney and gall bladder surgeries efficiently




③ Perfect X-Ray Imaging

- Head and leg sections with excellent radiolucent table guarantee high quality X-ray imaging
- Exchangeable head & leg sections
- Maximum 300 mm longitudinal shifting function

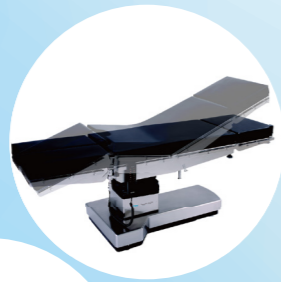


④ Surgical Configuration

- One-key reset function
- Flex and reflex function
- Flexible patient positionings with optional accessories are suitable for a wide variety of surgeries

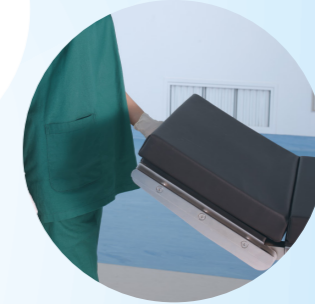

⑤ Ergonomic Adjustment

- Convenient leg plate and head plate positionings ensure one-hand control operation
- Imported gas springs guarantee smooth movement and accurate positioning
- Imported round-hole clamps enable easy lock and unlock accessories, ensuring safety during surgery time
- Wide range of table height adjustment and T-shape base provide comfortable working environment

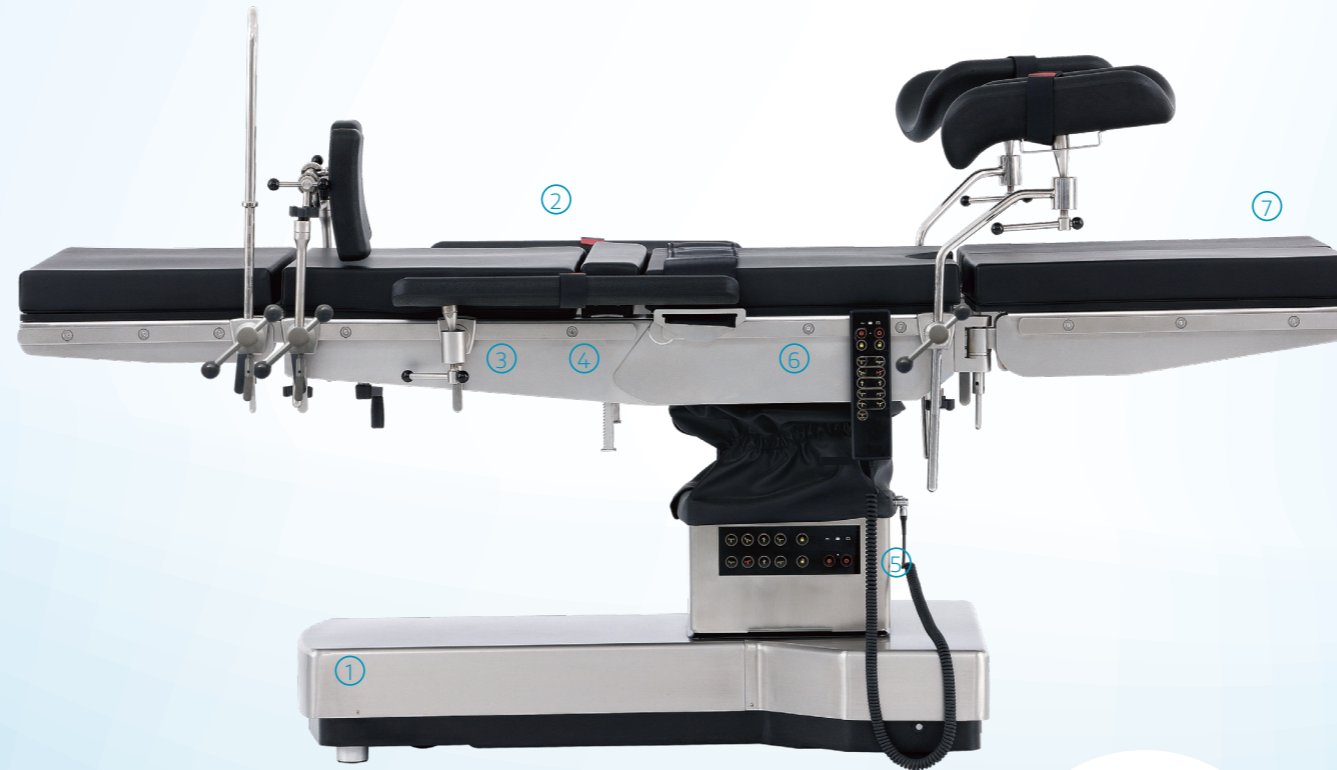

⑥ Safety

- Auto-lock system: table base will be locked automatically once the position was adjusted by controller, preventing possible misoperation
- Two control systems: manual and column control panel, allowing for backup if one is broken, ensuring the safety of surgery
- Intelligent safety system: Control panel will power off automatically when there is no operation over 1 minute, avoiding any possible misoperation



⑦ Pads

- High density memory foam with soft texture can effectively reduce the occurrence of bedsores













**Accessories**

- High strength 304 stainless steel
- Strong corrosion resistance, easy to clean and disinfection
- Simple structure, easy to install and disassemble
- Multiple options support multitude body positions and surgical procedures

**Standard Accessories**

		
Round-hole clamp	Anesthesia screen*	Body strap
		
Arm rest	Lateral support*	Knee crutch*

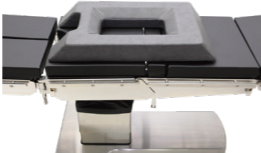**Optional Accessories**

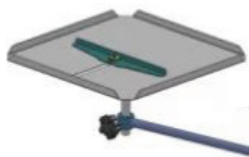		
Foot plate (notched clamp)**	High position arm holder*	Connection bracket
		
Head rest(connection bracket)***	Head plate(connection bracket)***	Hand operating table*
		
Instrument tray*	Iliac support*	Intervertebral pad
		
Drainage bowl	Infusion frame*	X-Ray Cassette Holder*



Orthopedic traction

**Orthopedic traction**




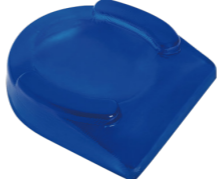
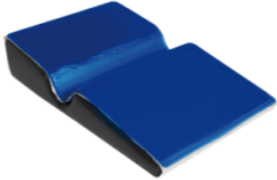




- In line with EU standards
- Suitable for most orthopedic surgeries
- High strength 304 stainless steel
- Structural stability design ensures maximum intraoperative safety
- X-ray radiolucent for free access to C-arm

**Noted:** \* means using together with the round-hole clamp  
 \*\* means using together with notched clamp  
 \*\*\* means using together with connection bracket

**Gel Pad**

All gel pads are made of viscoelastic polymer gel. They make patients feeling stable and comfortable during positioning and can be operated in the working environment of -10 C ~ 50 C .

- Uniform dispersion and biocompatibility of human skin and tissues to avoid pressure sores
- Easy to clean and disinfect
- X-ray translucent performance for X-ray imaging
- Cost effective and durable

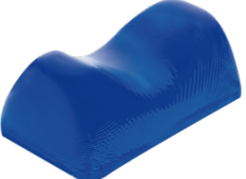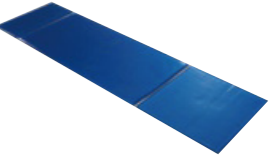		
Prostrate head pad	Donut head pad	Horseshoe head pad
		
Ophthalmic head pad	Lateral pad	Universal dome shaped pad
		
Contoured arm/Leg pad	Heel pad	Torso section

Flexible and Accurate Patient Positioning:



**Knee Arthroscopy Surgery**

Recommended accessories:  
 Donut head gel pad x1  
 Torso section gel pad x1  
 Arm rest x2  
 Knee support x1  
 Round-hole clamp x1



**Rectal Surgery**

Recommended accessories:  
 Arm rest x2  
 Foot plate x1 (pair)  
 Notched clamp x2  
 Iliac support x1  
 Round-hole clamp x2



**Chest/Abdominal Surgery**

Recommended accessories:  
 Donut head gel pad x1  
 Torso section gel pad x1  
 Foot plate x1 (pair)  
 Notched clamp x2  
 Arm rest x2



**E.N.T Surgery**

Recommended accessories:  
 Donut head gel pad x1  
 Torso section gel pad x1



**Orthopedic Surgery**

Recommended accessories:  
 Donut head gel pad x1  
 Torso section gel pad x1  
 Orthopedic traction x1  
 Arm rest x1



**Thoracic Kidney/ Gall Bladder Surgery**

Recommended accessories:  
 Donut head gel pad x1  
 Torso section gel pad x1  
 Arm rest x1  
 High position arm holder x1  
 Heel pad x1  
 Round-hole clamp x1



**Urology/OB/Gynecology Surgery**

Recommended accessories:  
 Donut head gel pad x1  
 Torso section gel pad x1  
 Knee crutch x2 Arm rest x2  
 Stirrup with gasspring adjustment for lithotomy surgery x1 (To protect peroneal nerve)  
 Drainage bowl x1



**Spinal Surgery**

Recommended accessories:  
 Prostrate head gel pad x1  
 Intervertebral pad x1  
 Arm rest x2  
 Torso section gel pad x1  
 Universal dome shaped pad x1



**Orthopedic Procedures**

Recommended accessories:  
 Donut head gel pad x1  
 Torso section gel pad x1  
 Orthopedic traction x1  
 Knee support for orthopedic traction x1  
 Installing rack of inserted rod x1  
 Arm rest x1



**Orthopedic Surgery for Lateral Femoral Extension**

Recommended accessories:  
 Orthopedic traction x1  
 Anti-tensile bracket for side-lying position x1  
 High position arm holder x1  
 Arm rest x1  
 Round-hole clamp x1



**Shoulder Surgery (Except Arthroscopic Surgery)**

Recommended accessories:  
 Arm rest x2  
 Head rest x1  
 Connection bracket x1  
 Body strap x1  
 Torso section gel pad x1



**Neurosurgery**

Recommended accessories:  
 Arm rest x2  
 Head rest x1  
 Connection bracket x1  
 Torso section gel pad x1  
 Universal dome shaped pad x1